



# WaveWare Technologies, Inc.

"We Deliver Information at the Speed of Light"



## Call Manager Software

### Call Manager Software: *Assign Window*

#### *For Wireless Call Devices*

Click on the "Assign" Button to View the Assign Wireless Call Devices Screen

Time	Log Item
3/12/13 2:14 PM	Loading Devices for Wired System Input Port...
3/12/13 2:14 PM	Loading Devices for Wireless System Input Port...
3/12/13 2:14 PM	Loading Stations...
3/12/13 2:14 PM	Loading Rooms/Residents...
3/12/13 2:14 PM	Loading Pagers...
3/12/13 2:14 PM	Loading Station/Pager Relationships...
3/12/13 2:14 PM	Loading Station/Room Relationships...
3/12/13 2:14 PM	Loading Device Types...
3/12/13 2:14 PM	Loading Receivers...
3/12/13 2:14 PM	Loading Report Templates...
3/12/13 2:14 PM	Logging Mode Activated
3/12/13 2:14 PM	Opened Log File - 031213.log
3/12/13 2:14 PM	Application Started - 031213.log
3/12/13 2:14 PM	Activating TAP on Input 1
3/12/13 2:14 PM	Opened COM Port 4, Settings: 9600,n,8,1 for Input 1
3/12/13 2:14 PM	Activating Inovonics on Input 2
3/12/13 2:14 PM	Opened COM Port 1, Settings: 9600,n,8,1 for Input 2
3/12/13 2:14 PM	Serial Port TAP Protocol Output Activated
3/12/13 2:14 PM	Opened COM Port 8, Settings: 9600,n,8,1 for TAP Protocol Serial Port Output
3/12/13 2:14 PM	System Successfully Started
3/12/13 2:14 PM	INITIALIZING INOVONICS (tm) RECEIVER ON COM 1
3/12/13 2:14 PM	Inovonics (tm) Receiver Initialized

No Activity on Input #2 - INOVONICS

This Feature is to allow a "User" the Ability to "Assign" and to "Un-Assign" or "Re-Assign" a Room to a Call Device, without the need of having Access to the Main Configuration Screens.

This Protects the System from becoming Corrupt by an Accidental Change to Important Settings/Configurations.

**NOTE:** You Must use the Password Feature in General Settings Step 14 (Page 37) for this to be Effective. You Do Not have to use Passwords, but it is Strongly Suggested that you do so.


Assign Window Continues on Next Page...

## Assign Window - Continued

### For Wireless Call Devices

In this Screen, you have the Ability to Change the Room Assignment.

You may “Un-Assign” or “Re-Assign” the Call Device to another Room Number / Resident.

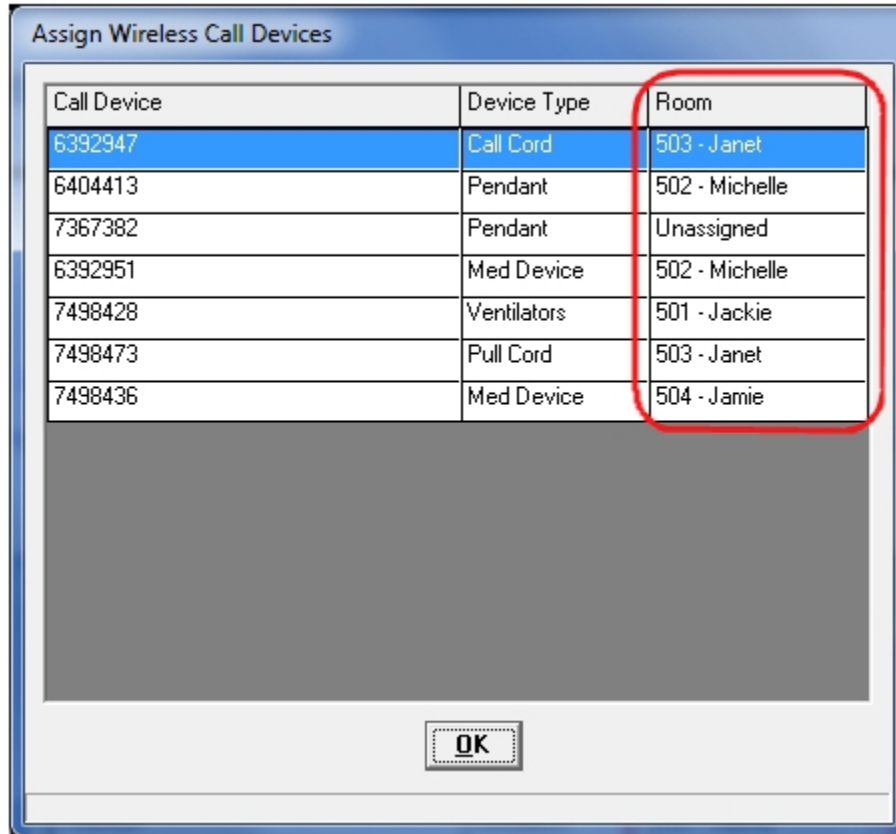


Call Device	Device Type	Room
6392947	Call Cord	503 - Janet
6404413	Pendant	502 - Michelle
7367382	Pendant	Unassigned
6392951	Med Device	502 - Michelle
7498428	Ventilators	501 - Jackie
7498473	Pull Cord	503 - Janet
7498436	Med Device	504 - Jamie

OK

Click on the “Room” Field and use the Drop down List to choose a “Room Number / Resident”.  
Rooms and Residents are Configured in Step 8 (Page 22).

If you leave a Room “Un-Assigned” you get this Warning Window. It is just a “Reminder”, you can Click “OK”.

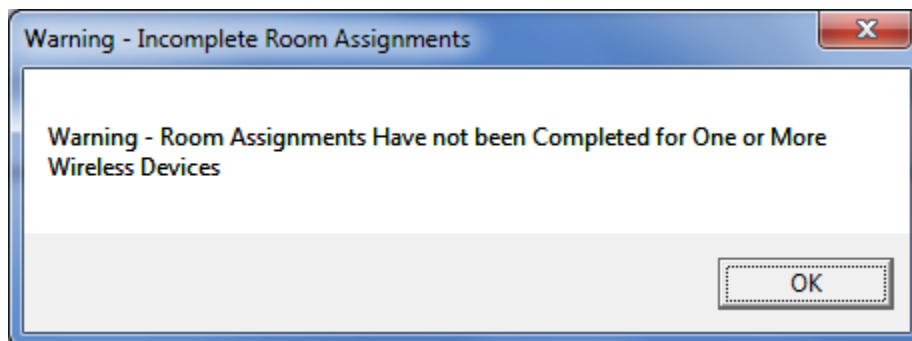

Warning - Incomplete Room Assignments

Warning - Room Assignments Have not been Completed for One or More Wireless Devices

OK

**NOTE:**

To Add a New Pendant or other Call Device, you would need to Use the Call Manager User Manual.  
See Step 9, of the Main Configuration Screens.

Please Contact your “Software Dealer” for Assistance or Questions.